


Plastic, Oral, Maxillofacial

April Product News

Cotswold Needles

Micro-Dissection Electrosurgery Needles

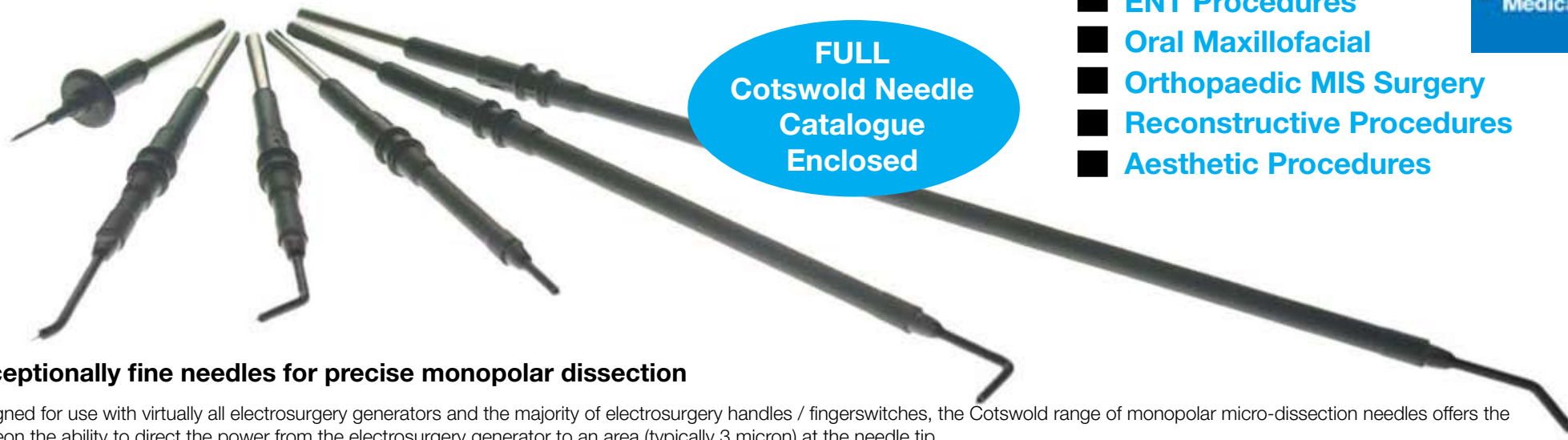

FULL
Cotswold Needle
Catalogue
Enclosed

Typical uses of the cotswold
needle electrodes are:

- ENT Procedures
- Oral Maxillofacial
- Orthopaedic MIS Surgery
- Reconstructive Procedures
- Aesthetic Procedures

Prima
Medical Limited

Exceptionally fine needles for precise monopolar dissection

Designed for use with virtually all electrosurgery generators and the majority of electrosurgery handles / fingerswitches, the Cotswold range of monopolar micro-dissection needles offers the surgeon the ability to direct the power from the electrosurgery generator to an area (typically 3 micron) at the needle tip.

With such a small and focussed point, the amount of power required to create a surgical effect is dramatically reduced as the electrosurgery current is concentrated to the smaller size – this reduction minimises the heat produced and hence the thermal energy spread.

- **Only low wattages of electrosurgery power required**
- **Reduced tissue necrosis**
- **Precision electrosurgery effect**
- **Reduced smoke**
- **Faster healing**

Safe - unique fingergrip allows for
easy insertion / extraction without
limiting vision.



Traditionally produced "sharp" electrosurgery needles become blunt very quickly during a procedure due to the stainless steel "melting", particularly with such a fine point. The Cotswold needles are produced in rigid, high temperature pure tungsten which has a melt temperature of 3400°C, significantly higher than those temperatures used in electrosurgery



Screw Implantation & Extraction Systems

Universal Mini Ratcheting Screwdriver System

Innovators of

Xtract-All

Judd Medical Limited is pleased to present the SHUKLA™ Medical range of innovative Orthopaedic Products.

The SHUKLA™ Medical products have been designed to offer a dedicated range of implantation & extraction systems offering Clinical Staff a wide variety of solutions for specific orthopaedic procedures. All systems are designed for **universal application** and have been uniquely developed to remove any manufacturers implants, giving greater flexibility in revision or trauma surgery.

Judd Medical Ltd
are **EXCLUSIVE**
UK DISTRIBUTORS
for SHUKLA™ Medical

DEMONSTRATIONS AVAILABLE !

Book your appointment with our
CUSTOMER SERVICE TEAM on
Tel +44(0)1527 559010



Installs & Removes:

- Small Screws typically encountered in:
 - Hand Surgery
 - Foot Surgery
 - Maxillofacial Surgery

Attachments:

- 24 Tips
 - Osteotome
 - 3 Torx® Drives
 - Phillips® Drive
 - Clutch Drive
 - 4 Centre Drives
 - 2 Slotted Drives
 - 3 Cruciform Drives
 - 3 Hex Drives
 - 4 Stripped Screw Extractors
 - 2 Broken Screw Extractors
 - Ratcheting Handle

Judd Medical Limited,

Highfield House, 53 Worcester Road, Bromsgrove, Worcestershire, B61 7DN, UK

Tel: +44 (0)1527 559010 Fax: +44 (0)1527 559080

Email: sales@judd-medical.co.uk

www.judd-medical.co.uk